

# WHEELCHAIRS


## BASIC MANUAL WHEELCHAIRS



## MARLIN

### Manual, foldable, steel wheelchair with quick release axles

Standard wheelchair on a cross frame made of chrome-molybdenum steel, equipped with quick release mechanism for rear wheels. It has height adjustable, lifting and detachable armrests. Its footrests are adjustable in length, rotating and removable. The wheelchair can also be height adjusted, thanks to special holes in the front fork and an assembly plate at the rear wheel. Seat height can be lowered from standard 50 cm to 48 cm, which is useful for shorter users or those after amputations. After replacing standard front wheels with smaller 6" size (optional equipment not included in standard price) the seat height can be lowered to 46 cm. Thanks to the rear wheel's assembly plate it is also possible to adjust seat angle. The Marlin can be attached easily to a caterpillar stairclimber.

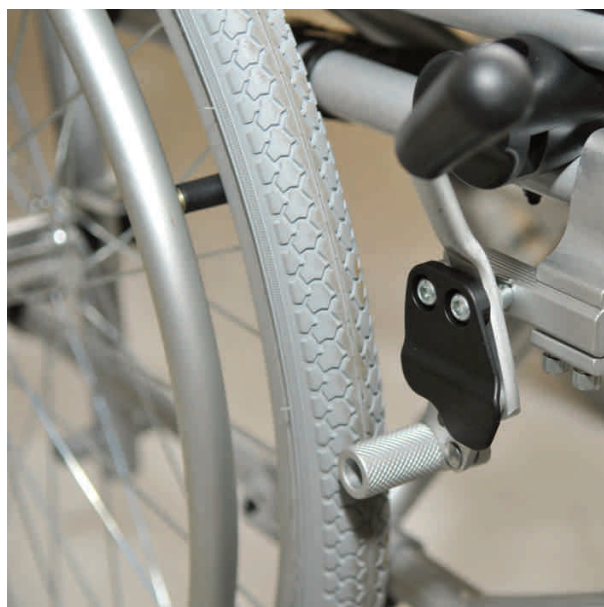

Front wheels size	Rear wheels size
200 x 37 mm (PU) – puncture proof	24" x 1 3/8" PU puncture proof tires Can be replaced with pneumatic tires

Seat width	Seat depth	Backrest height	Overall width	Seat to floor height	Seat to footplate height	Overall weight)	Max weight of user
40, 44, 48, 51 cm	41 cm	43 cm	61,5; 65,5; 69,5; 72,5; cm	50 cm	34-48 cm	21,6 kg	150 kg